

# **Quantum Analogs**

## **Chapter 2**

### **Advisor Manual**

# **Modeling a Hydrogen Atom**

## **with a**

### **Spherical Resonator**

**Professor Rene Matzdorf**  
**Universitaet Kassel**

## 2. Modeling a hydrogen atom with a spherical resonator

### Background:

The hydrogen atom, with a single electron in the Coulomb potential of the nucleus, is an ideal object for studying the basic principles of atomic physics. As the simplest of all atoms, without any electron correlations, it can be solved analytically.

The spherical symmetry of the three-dimensional problem makes it possible to separate the angular and radial variables for the solution of Schrödinger's equation. The acoustic analog uses a spherical resonator that allows a separation of variables for the solution of the Helmholtz equation in the same way as is done for the hydrogen atom. We will see that the eigenfunctions with respect to the angular variables – the spherical harmonics  $Y_l^m(\theta, \varphi)$  – are exactly the same for both problems. The radial eigenfunctions, however, are different.

The three-dimensional Schrödinger equation

$$E\psi(\vec{r}) = -\frac{\hbar^2}{2m}\Delta\psi(\vec{r}) - \frac{e^2}{r}\psi(\vec{r}) \quad (2.1)$$

expressed in polar coordinates

$$E\psi = \frac{\hbar^2}{2mr^2} \frac{\partial}{\partial r} \left( r^2 \frac{\partial \psi}{\partial r} \right) + \frac{\hbar^2}{2mr^2 \sin \theta} \frac{\partial}{\partial \theta} \left( \sin \theta \frac{\partial \psi}{\partial \theta} \right) + \frac{\hbar^2}{2mr^2 \sin^2 \theta} \frac{\partial^2 \psi}{\partial \varphi^2} - \frac{e^2}{r} \psi$$

can be separated in two differential equations with the ansatz

$$\psi(r, \theta, \varphi) = Y_l^m(\theta, \varphi) \chi_l(r). \quad (2.2)$$

The spherical harmonics are solutions of the differential equation

$$-\left[ \frac{1}{\sin \theta} \frac{\partial}{\partial \theta} \left( \sin \theta \frac{\partial}{\partial \theta} \right) + \frac{1}{\sin^2 \theta} \frac{\partial^2}{\partial \varphi^2} \right] Y_l^m(\theta, \varphi) = l(l+1) Y_l^m(\theta, \varphi) \quad (2.3)$$

and  $\chi_l(r)$  is a solution of the so called radial equation

$$-\frac{\hbar^2}{2mr} \frac{\partial^2}{\partial r^2} r \chi(r) - \frac{l(l+1)\hbar^2}{2mr^2} \chi(r) - \frac{e^2}{r} \chi(r) = E \chi(r). \quad (2.4)$$

In the case of the spherical acoustic resonator we transform eqn. 1.4

$$\frac{\partial^2 p}{\partial t^2} = \frac{1}{\rho \kappa} \Delta p \quad (2.5)$$

with the ansatz  $p(\vec{r}, t) = p(\vec{r}) \cos(\omega t)$  into the time independent Helmholtz equation

$$\omega^2 p(\vec{r}) = -\frac{1}{\rho \kappa} \Delta p(\vec{r}), \quad (2.6)$$

Using  $c$  as the speed of sound, equation 2.6 can be written as

$$-\frac{\omega^2}{c^2} p(\vec{r}) = \Delta p(\vec{r}) \quad (2.7)$$

The Helmholtz equation in polar coordinates is given by

$$-\frac{1}{r^2} \frac{\partial}{\partial r} \left( r^2 \frac{\partial p}{\partial r} \right) - \frac{1}{r^2 \sin \theta} \frac{\partial}{\partial \theta} \left( \sin \theta \frac{\partial p}{\partial \theta} \right) - \frac{1}{r^2 \sin^2 \theta} \frac{\partial^2 p}{\partial \varphi^2} = \frac{\omega^2}{c^2} p$$

It separates into a radial-function  $f(r)$  and the spherical harmonics  $Y_l^m(\theta, \varphi)$ .

$$p(r, \theta, \varphi) = Y_l^m(\theta, \varphi) f(r) \quad (2.8)$$

With this ansatz the Helmholtz equation is separated in one differential equation for the spherical harmonics

$$-\left[ \frac{1}{\sin \theta} \frac{\partial}{\partial \theta} \left( \sin \theta \frac{\partial}{\partial \theta} \right) + \frac{1}{\sin^2 \theta} \frac{\partial^2}{\partial \varphi^2} \right] Y_l^m(\theta, \varphi) = l(l+1) Y_l^m(\theta, \varphi) \quad (2.9)$$

and another for the radial function

$$-\frac{\partial^2 f}{\partial r^2} - \frac{2}{r} \frac{\partial f}{\partial r} + \frac{l(l+1)}{r^2} f(r) = \frac{\omega^2}{c^2} f(r) \quad (2.10)$$

You see immediately that eqn. 2.3 and eqn. 2.9 are exactly the same and have the same eigenfunctions and eigenvalues for the quantum numbers  $l$  (angular momentum or azimuthal quantum number) and  $m$  (magnetic quantum number). The radial equations are different, which, of course, results in different solutions. The Coulomb potential only appears in the radial equation (eqn. 2.4). Therefore, it does not affect the spherical harmonics. The eigenvalues of the radial equations are numerated by the quantum number  $n'$  (radial quantum number).

The energy levels  $E_{n'l}$  of the hydrogen atom are the eigenvalues of the radial equation (2.4) and the eigenfrequencies of the spherical acoustic resonator  $\omega_{n'l}$  are eigenvalues of the radial equation (2.10). Since the two differential equations are of different form, the resonance frequencies in the resonator can not be compared quantitatively with the energy levels of the hydrogen atom. However, the resonances can be classified with the same quantum numbers  $n'$  (radial quantum number),  $l$  (azimuthal quantum number) and  $m$  (magnetic quantum number). The quantum numbers are integers and

$$n' \geq 0 \quad l \geq 0 \quad -l \leq m \leq l \quad (2.11)$$

In the non-relativistic description of the hydrogen atom, many energy levels are degenerate, due to the special form of the Coulomb potential. The energies can be written in the form

$$E_{n'l} = - \left( \frac{e^2}{\hbar c} \right)^2 \frac{mc^2}{2(l+1+n')^2}. \quad (2.12)$$

All levels with the same value for  $(l+1+n')$  are degenerate. Therefore, a new quantum number is introduced that is called the “principal quantum number”  $n$ . It is given by


$$n = l + 1 + n' \quad (2.13)$$

For a given principal quantum number  $n$  the azimuthal quantum number  $l$  can take the values

$$0 \leq l \leq n-1 \quad (2.14)$$

even though it runs to infinity for a given radial quantum number.

In the diagrams of the hydrogen atom spectrum shown below, the energy levels are labeled in two different ways. In the left figure they are labeled in the ordinary manner, using the principal quantum number. The right figure shows the energy levels labeled using the radial quantum number.

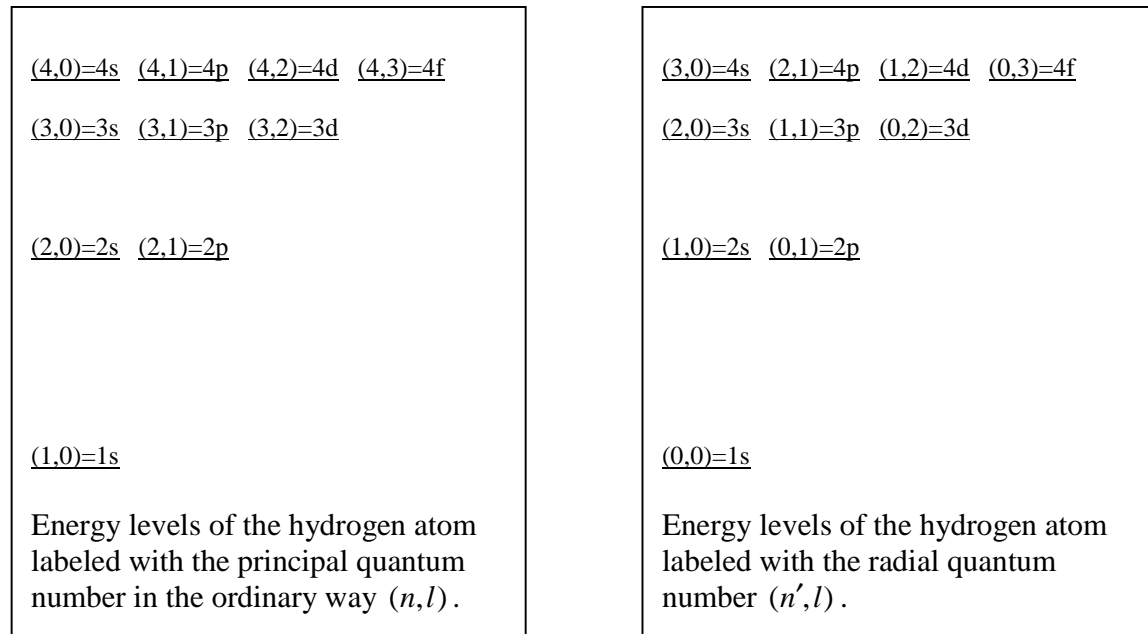

The degeneracy of levels with the same principal quantum number does not have an analog in the spherical acoustic resonator, since the radial equation is different.

In the spherically symmetric case, the eigenvalues for different magnetic quantum numbers  $m$  are degenerate for any form of the radial equation. This is true for both the hydrogen atom and the spherical acoustic resonator. In general, the eigenvalues numbered by the quantum numbers  $(n,l)$  or by  $(n',l)$  are  $(2l+1)$ -fold degenerate. This degeneracy is lifted when the spherical symmetry is broken.

Now let's do some experiments that allow us to see many of these effects. First, we will identify the resonances by their angular dependence.

## 2.1 Measure resonances in the spherical resonator and determine their quantum numbers

**Objective:** Determine the resonance frequencies for the spherical resonator and gather data to determine their angular-momentum quantum numbers.

### Equipment Required:

TeachSpin Quantum Analog System: Controller, Hemispheres, Accessories


Sine wave generator capable of producing 1-50 kHz with a peak-to-peak voltage of 0.50 V

Two-Channel Oscilloscope

### Setup:

Assemble two of the hemispheres so that the speaker is in the lower hemisphere and a microphone is in the upper hemisphere. (Looking carefully at the photo, you will see the speaker wire at the lower right.) Adjust the position of the upper hemisphere so that  $\alpha = 180^\circ$  on the scale is at the reference mark. In this position, the speaker and the upper microphone are at opposite ends of a diameter. (The microphones will be one above the other.)

Attach a BNC splitter or “tee” to *SINE WAVE INPUT* on the Controller. Connect the output of your sine wave generator to one side of the splitter. Use a BNC cable to send the sound signal to Channel 1 of the oscilloscope. Plug the lead from the speaker on the lower hemisphere to *SPEAKER OUTPUT* on the Controller. The same sine wave now goes to both the speaker and Channel 1.

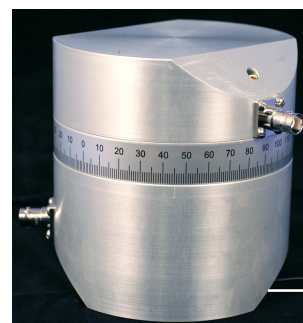

**Atom Analog**  
Microphones are imbedded  
under BNC connectors.  
Speaker is at lower right.

Use a BNC cable to connect the microphone output from the **upper** hemisphere to *MICROPHONE INPUT* on the Controller. Connect *AC MONITOR* on the Controller to Channel 2 of the oscilloscope to display the sound signal received by the microphone. Trigger the oscilloscope on Channel 1.

Use the *ATTENUATOR* dial on the Controller to keep the signal on Channel 2 from going off scale. Remember, with an attenuator, a higher reading on the dial gives a smaller signal. (Appendix 1 describes the function of each part of the Controller.)

### Experiment:

Start at a low frequency and sweep the frequency up to about 8 kHz (8000 Hz).

Write down all the resonance frequencies you observe. (If you listen carefully, you may actually hear some of them.)

### ADVISOR INFORMATION:

At this angle between speaker and microphone, the resonances with high angular quantum number  $l$  are visible with intensities higher than those observed at any other angle. Depending on the temperature, you will find resonances near the following frequencies: 2291 Hz, 3678 Hz, 4968 Hz, 6210 Hz, 6539 Hz, 7425 Hz, 8023 Hz ...

**Objective:**

Observe, qualitatively, the way the amplitude of the resonance signal depends on the location of the microphone.

**Experiment:**

We will now gather data that will allow us to infer the angular quantum numbers of the resonances. Go to the second resonance, at about 3680 Hz. Fine-tune the frequency until it is as close as possible to the peak of the resonance. Shift the curves on the oscilloscope horizontally so that a maximum of the microphone signal (Channel 2) is located in the center of the image and marked by a vertical line. Now, watching the signal on the oscilloscope, slowly rotate the upper hemisphere, with respect to the lower one, from  $\alpha = 180^\circ$  to  $\alpha = 0^\circ$ .

**Questions:**

How did the amplitude change? Did the signal change sign? Find the angle where the amplitude is zero. At which angles is the signal largest? Do both extrema have the same amplitude?

**Note:** Do not warm the aluminum parts too much by touching them with your hands. The speed of sound is temperature-dependent, and, in consequence, the resonance frequency would shift with temperature. While analyzing the angular dependence, the chosen generator frequency should remain on top of the resonance.


**ADVISOR INFORMATION:**

The amplitude changes sign and has a node at about  $\alpha = 100^\circ$  corresponding to the node in the Legendre Polynomial at  $\theta = 125.26^\circ$ . The maximum at  $\alpha = 0^\circ$  ( $\theta = 90^\circ$ ) has half the magnitude of the maximum at  $\alpha = 180^\circ$  ( $\theta = 180^\circ$ ). Experimentally, the results are quite sensitive to the adjustment of frequency and the positioning of the upper hemisphere. Pressing gently on the upper hemisphere will give you a feeling for the sensitivity in the measurement of these amplitudes.


**Analyze the data:**

The angle  $\alpha$  read on the scale is not a suitable angle for comparison with theory. Notice that the scale reading,  $\alpha$ , running from  $0 - 180^\circ$ , tells you the rotation of the upper hemisphere about a vertical axis. The symmetry axis for this system, however, is determined by the speaker. The angle of interest, therefore, is measured using the speaker location as zero. To analyze the data, you must first use  $\alpha$  to calculate the polar angle,  $\theta$ . It is this angle which we use for polar coordinates. To clarify how this works try the following.

Assemble the sphere with the upper hemisphere set so that  $\alpha = 180^\circ$ . Temporarily open the resonator and notice that at this setting the speaker and the upper microphone are 180 degrees apart in space;  $\theta = 180^\circ$ . Reassemble the spheres, turn the upper sphere to  $\alpha = 0^\circ$ , and open the resonator again. You will see that the spatial separation of the speaker and microphone, the polar angle of interest, is now  $\theta = 90^\circ$ .

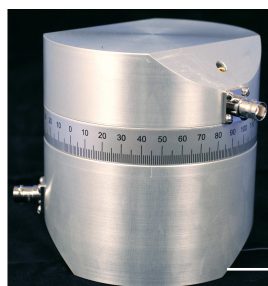

**Atom Analog**  
Speaker is at lower right

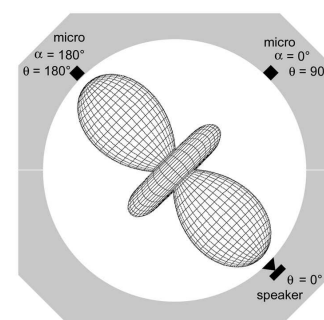

**Sample Sound Amplitude Pattern**

Both the speaker and microphone are at an angle of  $45^\circ$  with respect to the horizontal plane between the hemispheres. By rotating the hemispheres with respect to each other, the angle  $\theta$  can be changed from  $\theta = 90^\circ$  (at  $\alpha = 0^\circ$ ) to  $\theta = 180^\circ$  (at  $\alpha = 180^\circ$ ). Intermediate angles can be calculated using the formula

$$\theta = \arccos\left(\frac{1}{2} \cos \alpha - \frac{1}{2}\right). \quad (2.15)$$

You have measured the  $\theta$ -dependence of the spherical harmonic function  $Y_l^m(\theta, \varphi)$  with  $l = 2$  and  $m = 0$ . Now we need to learn more about the spherical harmonics to compare the experiment with theory.

#### Derivation of equation 2.15

Assume that the speaker is located in the x-z-plane and the vertical axis is the z-axis. The position of the speaker in a sphere with unit-radius is given by the vector  $\vec{s} = (\sqrt{\frac{1}{2}}, 0, -\sqrt{\frac{1}{2}})$ .

We want to calculate the angle between speaker and microphone, which is the angle  $\theta$ .

To calculate  $\theta$  we use rotary matrices. In the first step, we rotate the vector  $\vec{s}$  from the position of the speaker (vector  $\vec{s}$ ) by 90 degrees around the y-axis arriving at  $(\sqrt{\frac{1}{2}}, 0, \sqrt{\frac{1}{2}})$ .

In the second step, we rotate by the angle  $\alpha$  around the z axis. Lets call the resulting vector,  $\vec{m}$ , the position of the microphone.

From the scalar-product,  $\vec{m} \cdot \vec{s} = |\vec{m}| |\vec{s}| \cos \theta = \cos \theta$ , we get the angle  $\theta$ .

First rotation:

$$\begin{bmatrix} \cos 90^\circ & 0 & -\sin 90^\circ \\ 0 & 1 & 0 \\ \sin 90^\circ & 0 & \cos 90^\circ \end{bmatrix} \cdot \begin{bmatrix} \sqrt{\frac{1}{2}} \\ 0 \\ -\sqrt{\frac{1}{2}} \end{bmatrix} = \begin{bmatrix} \sqrt{\frac{1}{2}} \\ 0 \\ \sqrt{\frac{1}{2}} \end{bmatrix}$$

Second rotation:

$$\begin{bmatrix} \cos \alpha & -\sin \alpha & 0 \\ \sin \alpha & \cos \alpha & 0 \\ 0 & 0 & 1 \end{bmatrix} \cdot \begin{bmatrix} \sqrt{\frac{1}{2}} \\ 0 \\ \sqrt{\frac{1}{2}} \end{bmatrix} = \begin{bmatrix} \sqrt{\frac{1}{2}} \cos \alpha \\ \sqrt{\frac{1}{2}} \sin \alpha \\ \sqrt{\frac{1}{2}} \end{bmatrix}$$

scalar-product:


$$\vec{m} \cdot \vec{s} = \begin{bmatrix} \sqrt{\frac{1}{2}} \cos \alpha \\ \sqrt{\frac{1}{2}} \sin \alpha \\ \sqrt{\frac{1}{2}} \end{bmatrix} \cdot \begin{bmatrix} \sqrt{\frac{1}{2}} \\ 0 \\ -\sqrt{\frac{1}{2}} \end{bmatrix} = \frac{1}{2} \cos \alpha - \frac{1}{2}$$

with

$$\vec{m} \cdot \vec{s} = \cos \theta$$

we get the result:

$$\theta = \arccos\left(\frac{1}{2} \cos \alpha - \frac{1}{2}\right)$$

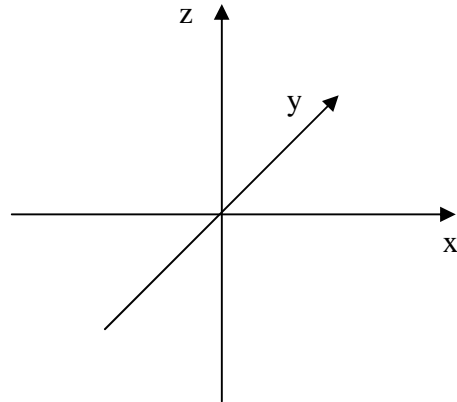

**Spherical Harmonics and Legendre Polynomials:**

The spherical harmonics  $Y_l^m(\theta, \varphi)$  can be written as

$$Y_l^m(\theta, \varphi) \propto P_l^m(\cos \theta) e^{im\varphi} \quad (2.16)$$

in terms of the associated Legendre polynomials  $P_l^m$ . For these experiments, we can restrict ourselves to the case  $m = 0$ , because our speaker creates waves with cylindrical symmetry about the speaker axis. For  $m = 0$  the spherical harmonics do not have a  $\varphi$ -dependence and the wave function has the same amplitude for all azimuthal angles,  $\varphi$ . The dependence on the polar angle  $\theta$  is given by the Legendre polynomials

$$Y_l^0(\theta, \varphi) \propto P_l^0(\cos \theta) \quad (2.17)$$

The first nine Legendre polynomials are shown below:

$$P_0(\cos \theta) = 1$$

$$P_1(\cos \theta) = \cos \theta$$

$$P_2(\cos \theta) = \frac{1}{2}(3 \cos^2 \theta - 1)$$

$$P_3(\cos \theta) = \frac{1}{2}(5 \cos^3 \theta - 3 \cos \theta)$$

$$P_4(\cos \theta) = \frac{1}{8}(35 \cos^4 \theta - 30 \cos^2 \theta + 3)$$


$$P_5(\cos \theta) = \frac{1}{8}(63 \cos^5 \theta - 70 \cos^3 \theta + 15 \cos \theta)$$

$$P_6(\cos \theta) = \frac{1}{16}(231 \cos^6 \theta - 315 \cos^4 \theta + 105 \cos^2 \theta - 5)$$

$$P_7(\cos \theta) = \frac{1}{16}(429 \cos^7 \theta - 693 \cos^5 \theta + 315 \cos^3 \theta - 35 \cos \theta)$$

$$P_8(\cos \theta) = \frac{1}{128}(6435 \cos^8 \theta - 12012 \cos^6 \theta + 6930 \cos^4 \theta - 1260 \cos^2 \theta + 35)$$

In Fig. 2.1 and 2.2 the first six Legendre polynomials are plotted. The number of nodes in each Legendre polynomial is equal to the azimuthal quantum number  $l$ .

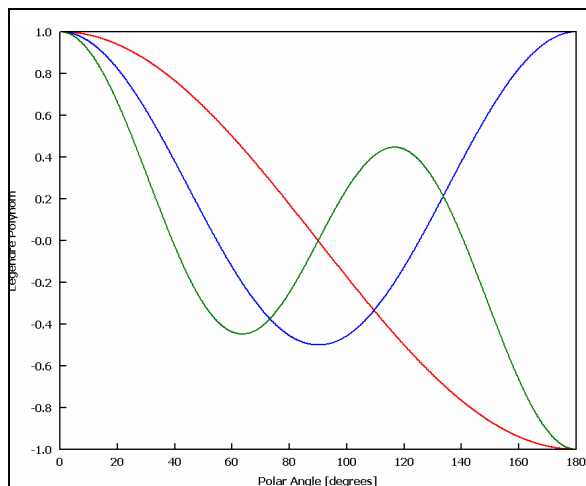


**Fig. 2.1:** Legendre Polynomials

$$P_1(\cos \theta) = \cos \theta \text{ in red,}$$

$$P_2(\cos \theta) = \frac{1}{2}(3\cos^2 \theta - 1) \text{ in blue and}$$

$$P_3(\cos \theta) = \frac{1}{2}(5\cos^3 \theta - 3\cos \theta) \text{ in green.}$$

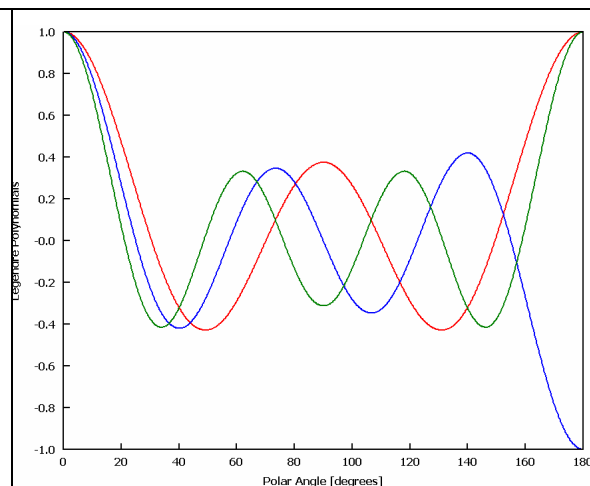

**Fig. 2.2:** Legendre Polynomials

$$P_4(\cos \theta) = \frac{1}{8}(35\cos^4 \theta - 30\cos^2 \theta + 3) \text{ in red,}$$

$$P_5(\cos \theta) = \frac{1}{8}(63\cos^5 \theta - 70\cos^3 \theta + 15\cos \theta) \text{ in blue}$$

$$P_6(\cos \theta) = \frac{1}{16}(231\cos^6 \theta - 315\cos^4 \theta + 105\cos^2 \theta - 5) \text{ in green.}$$

In the following table, the nodes of the Legendre polynomials are listed. **Be aware that these are the polar angles  $\theta$ , and not the angles you read on the scale.**

$P_0$								
$P_1$	$90^\circ$							
$P_2$	$54.74^\circ$	$125.26^\circ$						
$P_3$	$39.23^\circ$	$90^\circ$	$140.77^\circ$					
$P_4$	$30.56^\circ$	$70.12^\circ$	$109.88^\circ$	$149.44^\circ$				
$P_5$	$25.02^\circ$	$57.42^\circ$	$90^\circ$	$122.58^\circ$	$154.98^\circ$			
$P_6$	$21.18^\circ$	$48.61^\circ$	$76.19^\circ$	$103.81^\circ$	$131.39^\circ$	$158.82^\circ$		
$P_7$	$18.36^\circ$	$42.14^\circ$	$66.06^\circ$	$90^\circ$	$113.94^\circ$	$137.86^\circ$	$161.64^\circ$	
$P_8$	$16.20^\circ$	$37.19^\circ$	$58.30^\circ$	$79.43^\circ$	$100.57^\circ$	$121.70^\circ$	$142.81^\circ$	$163.80^\circ$

**Table 2.1:** Nodes of the first eight Legendre polynomials, given in the polar angles  $\theta$ .

**Questions:**

Now you can identify the angular quantum number  $l$  of the second resonance you have measured.

As you varied  $\alpha$  from  $180^\circ$  to  $0^\circ$ , what range of  $\theta$  did you cover?

How many nodes did you discover in the range you covered?

Based on your observations, to what  $l$  value does the resonance you examined correspond?

Does the  $\theta$  angle measurement of the node you have measured agree with the angle predicted by the theory?

Do the relative magnitudes of the extrema fit to the theory?

**Note about the magnetic quantum number:**

The resonance that you have analyzed is  $(2l+1)$ -fold degenerate with respect to the magnetic quantum number  $m$ . However, in this experiment we observe almost exclusively the  $m = 0$  state. The standing sound wave in the sphere is driven by the local speaker. The speaker defines the  $z$ -axis of the problem. It emits a wave traveling more or less back and forth along the  $z$ -axis and having cylindrical symmetry around that axis. This symmetry of the standing wave is described by the  $m = 0$  state. States with other  $m \neq 0$  describe waves that move on an orbit inside the sphere. These types of waves are much less effectively driven by our speaker located on the  $z$ -axis, since these states have nodes at  $\theta = 0^\circ$  and  $\theta = 180^\circ$ .

**Objective:** We will trace out the angular dependence of the amplitude of the wave function.

**Additional Apparatus:** dc voltmeter

**Setup:**

As in the first part of this experiment, attach a BNC splitter to *SINE WAVE INPUT* on the Controller. Connect the output of your sine wave generator to one side of the splitter. Use a BNC cable to send the sound signal to Channel 1 of the oscilloscope. Plug the lead from the speaker on the lower hemisphere to *SPEAKER OUTPUT* on the Controller. The same sine wave now goes to both the speaker and Channel 1.

Use a BNC cable to connect the microphone output from the upper hemisphere to *MICROPHONE INPUT*. Connect *AC MONITOR* on the Controller to Channel 2 of the oscilloscope to display the sound signal received by the microphone. Trigger the oscilloscope on Channel 1.

This time put the upper hemisphere in the position  $\alpha = 0^\circ$  on the scale. In this position the microphone is directly above the speaker which means angle  $\theta$  will be  $90^\circ$ .

To observe the amplitude of the sound signal at the microphone, connect a voltmeter to *DETECTOR-OUTPUT*. You should also observe the sound signal itself by connecting the *AC-MONITOR* on the Controller with Channel 2 of the oscilloscope. Use the *ATTENUATOR* dial on the Controller to keep the signal from going off scale. Trigger the oscilloscope to Channel 1.

**Experiment:**

For a couple of major resonances, measure the amplitude as function of the angle  $\alpha$ . You can read the absolute value of the amplitude on the voltmeter and use the oscilloscope to determine the sign.

Record the nodes (angle at which the amplitude is zero) for the same resonances.

**Analyze the data:**

Plot your data as function of the polar angle  $\theta$  and fit the data with the Legendre polynomial that is the best match. Do this for all the resonances you have measured.


Compare the nodes you have measured with the nodes of the corresponding Legendre polynomial given in table 2.1.

**Note:**


Some of the resonances are very close to each other so that the peaks are overlapping. This will result in a superposition of two wavefunctions with different quantum numbers. In this case, the angular dependence you have measured does not fit to a single Legendre polynomial. We will analyze these cases in more detail by taking spectra with the computer.

**ADVISOR INFORMATION:** These plots are for the advisor only.

Angular dependence of the amplitude for the resonance at 3678 Hz (open circles) and fit with the Legendre Polynomial  $P_2(\cos\theta)$  (solid line).

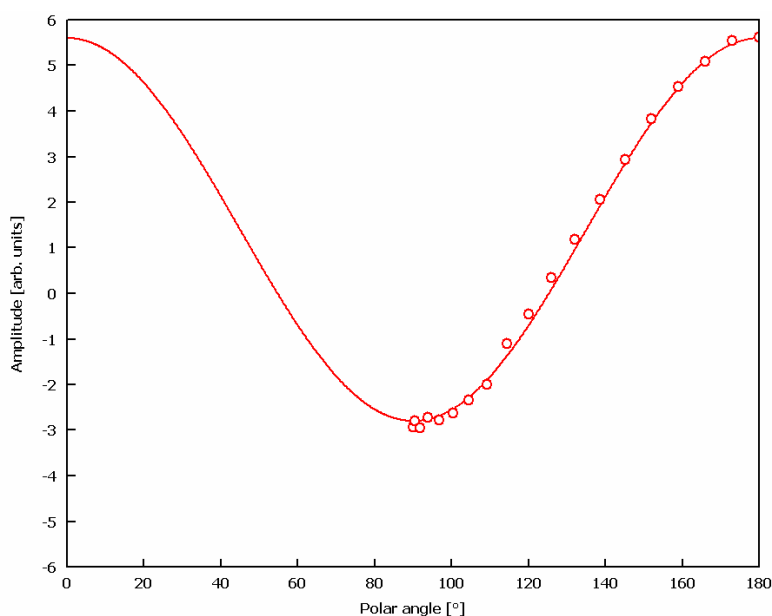


Angular dependence of the amplitude for the resonance at 4968 Hz (open circles) and fit with the Legendre Polynomial  $P_3(\cos\theta)$  (solid line).

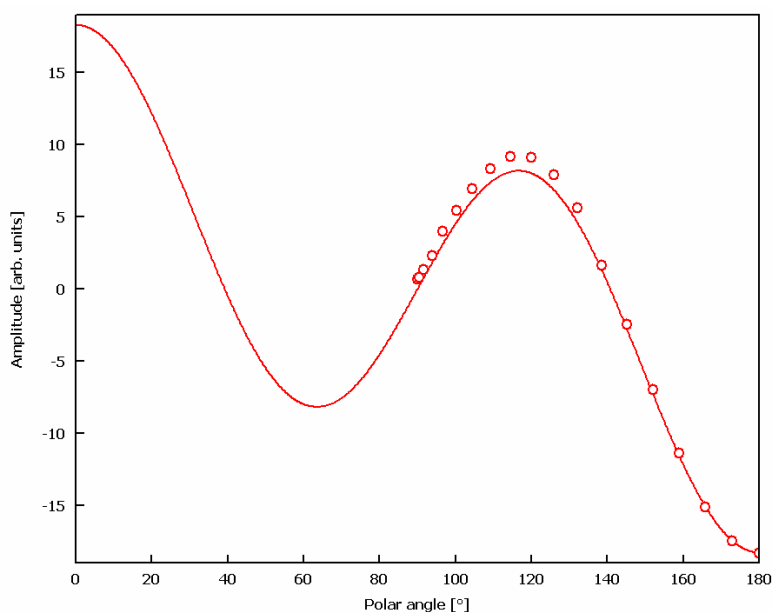


## ADVISOR INFORMATION CONTINUED

Angular dependence of the amplitude for the resonance at 7425Hz (open circles) and fit with the Legendre Polynomial  $P_5(\cos\theta)$  (solid line).

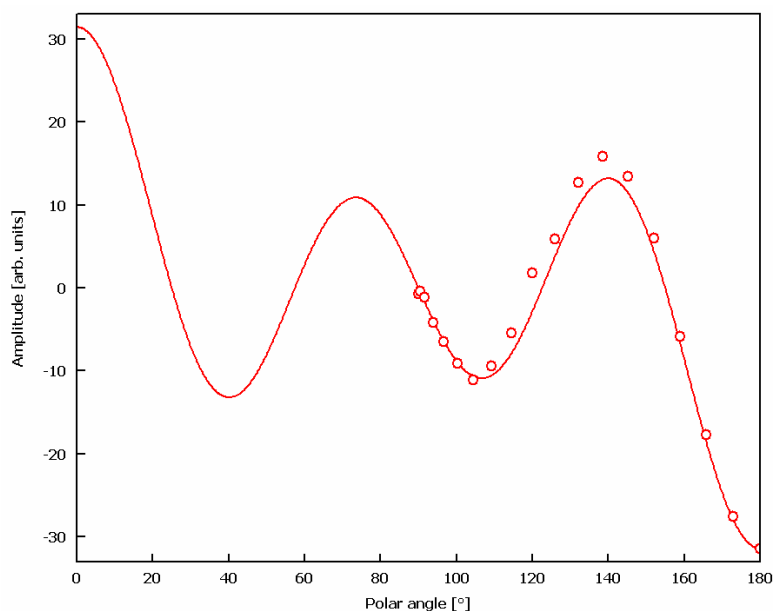

Angular dependence of the amplitude for the resonance at 9820Hz (open circles) and fit with the Legendre Polynomial  $P_7(\cos\theta)$  (solid line).

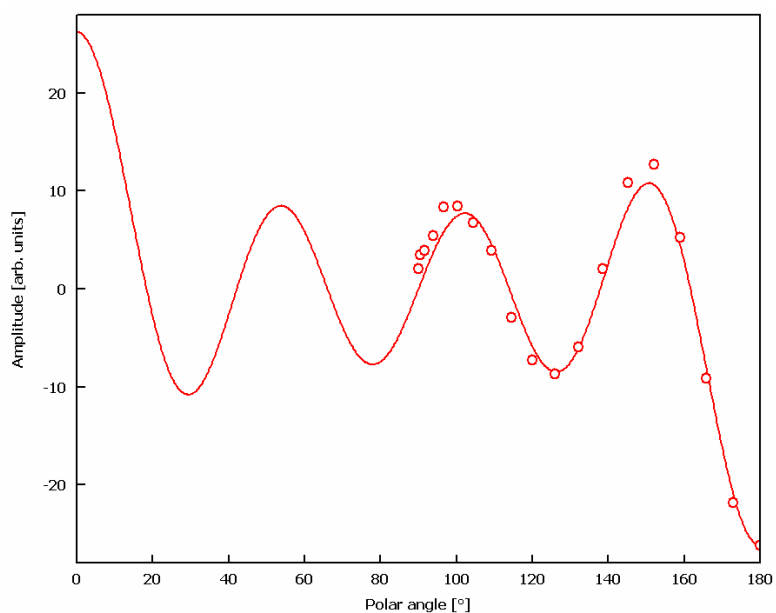

## 2.2 Measure spectra and wavefunctions in the spherical resonator with the computer

**Objective:** In this experiment, you will use a computer sound card both to generate the sound wave and to sweep its frequency. You will use the oscilloscope to observe the actual sine wave signals both going into the speaker and coming from the microphone. Simultaneously, you will use the computer to display a spectrum which shows the amplitude of the signal from the microphone as a function of the frequency of the sound.

### Equipment Required:

TeachSpin Quantum Analog System: Controller, Hemispheres, Accessories

Two-Channel Oscilloscope

Two adapter cables (BNC - 3.5 mm plug)

Computer with sound card installed and Quantum Analogs "SpectrumSLC.exe" running

**WARNING:** The BNC-to-3.5-mm adapter cables are provided as a convenient way to couple signals between the Controller and sound card. Unfortunately, they could also provide a way for excessive external voltage sources to damage a sound card. Most sound cards are somewhat protected against excessive inputs, but *it is the user's responsibility to ensure that adapter cable voltages are kept BELOW 5 Volts peak-to-peak.*

The maximum peak-to-peak value for optimum performance of the Quantum Analogs system depends on your sound card and can vary from 500 mV to 2 V.

### Setup:

Now, using connectors on the Controller, you will send the sound card signal to both the speaker and Channel 1 of the oscilloscope, and the microphone signal to both the microphone input of the computer and to Channel 2 of the oscilloscope.

**First, make sure that the ATTENUATOR knob on the Controller is set at 10.0 (out of 10) turns.**

Let's start with the sound signal. Attach a BNC splitter of "tee" to *SINE WAVE INPUT* on the Controller. Using the adapter cable, connect the output of the sound card to one arm of the splitter. With a BNC cable, convey the sound card signal from the splitter to Channel 1 of your oscilloscope. Plug the lead from the speaker on the lower hemisphere to *SPEAKER OUTPUT* on the Controller. The sound card signal is now going to both the speaker and Channel 1.

The microphone signal will also be sent two different places. Connect the microphone on the upper hemisphere to *MICROPHONE INPUT* on the Controller. Put a BNC splitter on the Controller connector labeled *AC-MONITOR*. From the splitter, use an adapter cable to send the microphone signal to the microphone input on the computer sound card. Use a BNC cable to send the same signal to Channel 2 of the oscilloscope to show the actual signal coming from the microphone.

The computer will plot the instantaneous frequency generated by the sound card on the x-axis and the amplitude of the microphone input signal on the y-axis.

**The next job is to adjust the magnitude of both the speaker and microphone signals so that you will have maximum signal while keeping the microphone input to the computer from saturating.** Peak-to-peak signals to the microphone input can range from 0.50 to 2.0 volts depending upon your sound card.

Once the program, SpectrumSLC.exe., is running, you can configure the computer. Go to the menu at the top of the screen and choose Configure > Input Channel/Volume. At this point, choose *Line In*, if it is available; otherwise choose *Microphone*. On this screen, set the microphone volume slider to the middle of its range.

To set the speaker volume, use the *Amplitude Output Signal* on the lower left of the computer screen. That slider should also be set to middle range.

The microphone signal coming from the apparatus first passes through a built-in amplifier, and then through the *ATTENUATOR*, before reaching the *AC-MONITOR* connector. The ten-turn knob on the attenuator *decreases* the incoming signal by a factor ranging from zero to 100. For example, a setting of 9.0 turns (out of the 10 turns possible) stands for an attenuation of 9/10 or 90% attenuation of the signal. (A higher setting means a smaller signal.)

After taking an initial wide range spectrum, choose a section that includes the highest peak and a smaller one next to it. Readjust the scan to cover just this portion. Using the option that allows you to keep successive spectra visible, take Spectrum 1, 2, 3, etc. with the attenuator knob set at 9.9, 9.8, 9.7, . . . turns (out of ten). The nesting heights of the peaks will tell you whether or not the system is behaving in a linear fashion. Continue to go lower on the 10-turn dial setting until the computer program flashes 'saturation'. You will also have visual evidence of saturation – a flat section on the tallest peak or a smaller “nesting” spacing. (See Appendix 2 or 3 for details.)

Once you have reached saturation, drop back into the linear range. Now you can operate with confidence that the signals you see really are proportional to the amplitude of the sound wave you are studying.

### **Experiment:**

Set the hemispheres so that the scale angle  $\alpha = 180^\circ$ .

Start the program SpectrumSLC.exe and measure an overview spectrum. You can use coarse steps such as 10 Hz and a short time per step such as 50 ms.


Change the angle between the upper and the lower hemisphere several times and observe the how the spectrum changes. Be sure to look at the spectrum for  $\alpha = 0^\circ$ .

**Question:** What changes do you notice?


**ADVISOR INFORMATION:**

The amplitudes of the peaks depend strongly on the angle between speaker and microphone due to the different symmetry of the wave functions. With the hemispheres set so that the scale angle  $\alpha = 180^\circ$  (which corresponds to  $\theta = 180^\circ$ ), the resonances with high  $l$  are at their most intense. Some of the resonances are only visible when another resonance nearby is suppressed by measuring in a position where its wave function has a node.

Spectrum measured in the  $\alpha = 0^\circ$  position of the hemispheres giving  $\theta = 90^\circ$  between speaker and microphone.

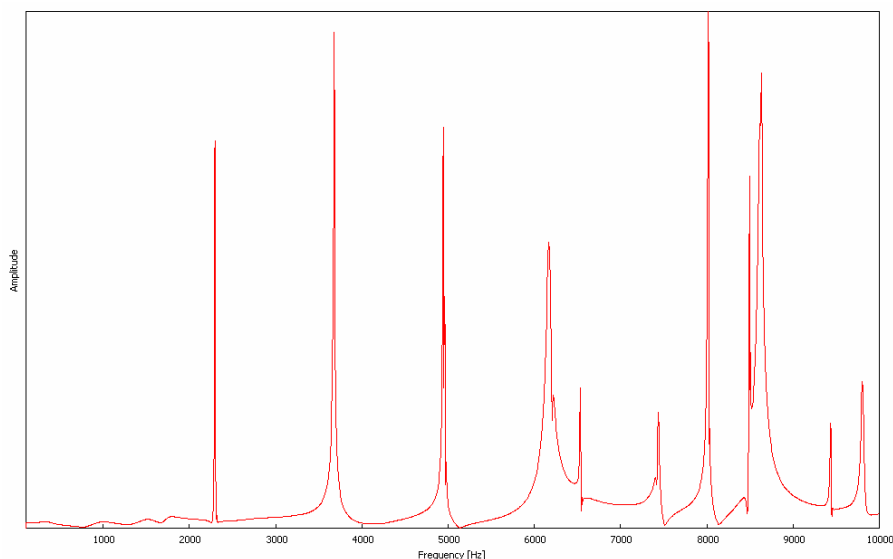

Spectrum measured in the  $\alpha = 180^\circ$  position of the hemispheres, giving  $\theta = 180^\circ$  between speaker and microphone.

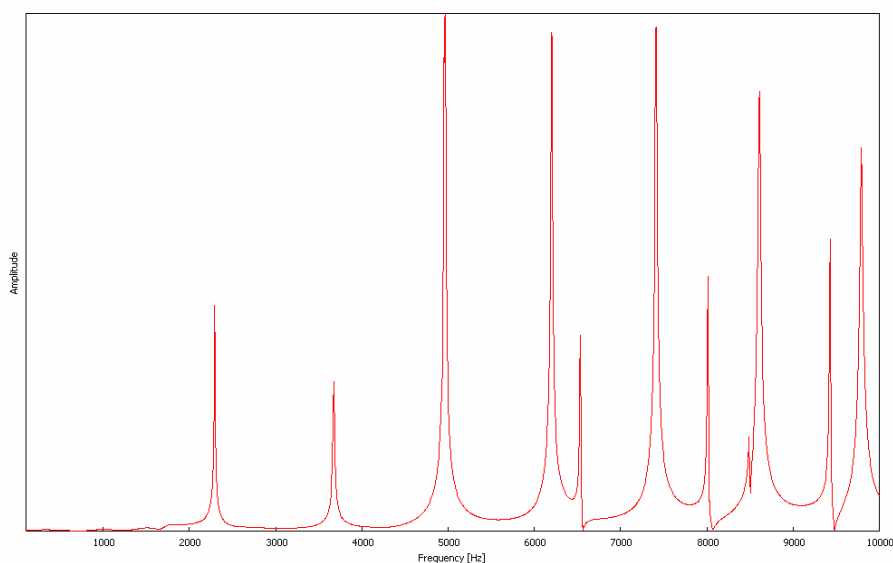


**Experiment:**

Go back to  $\alpha = 0^\circ$  and look in more detail at the peak near 5000 Hz. Actually, there are two peaks close to each other. Take a spectrum that measures slow enough and with sufficiently small steps to show the details of these two peaks. Also, take spectra for this range at  $\alpha = 20^\circ$  and  $\alpha = 40^\circ$ .

**Question:** What do you notice?

**ADVISOR INFORMATION:**

The double peak structure at about 5000 Hz belongs to two different resonances with different “angular momentum” values. The left peak is the resonance  $n' = 1, l = 0$ , the right one is the resonance  $n' = 0, l = 3$ . The amplitudes of the peaks depend differently on the angle  $\theta$  due to their different  $l$  values. For  $\theta = 90^\circ$ , the right peak has a minimum (node in the wave function) so that the left one is at its most visible.

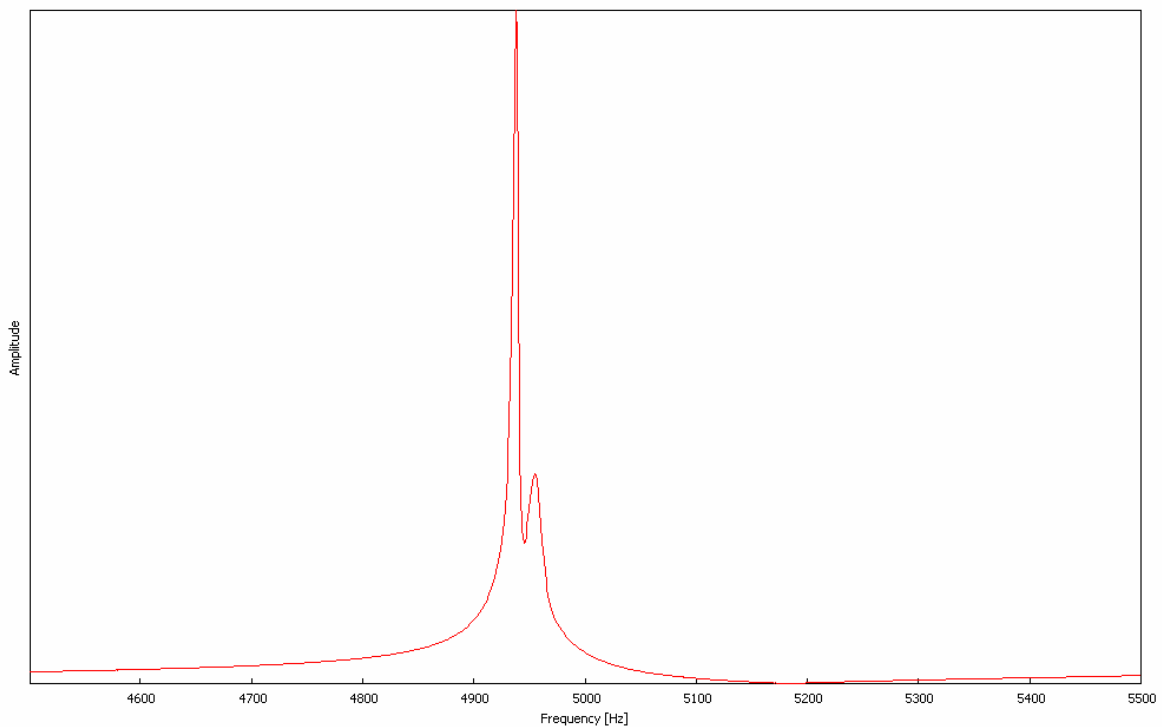

Spectrum of peaks near 5000 Hz measured with 1 Hz steps and 100 ms per step.

**Objective:** Create polar plots for a series of resonances and use the plots to identify the angular momentum number and spherical harmonic function of each resonance.

**Experiment:**

Now we will measure the wavefunctions of the different resonances and visualize them by a polar plot of the amplitude  $A(\theta)$ . The computer calculates the polar angle  $\theta$  from the angle  $\alpha$  and it plots the absolute value of the amplitude as function of  $\theta$  in a polar plot. This diagram makes it easy to identify the angular quantum number and the spherical harmonic function.

Take a spectrum with  $\alpha = 180^\circ$  from 2000 Hz to 7000 Hz sufficiently slowly. If you click with the left mouse button on a peak, the output frequency is adjusted to the value at which you clicked. Look at the oscilloscope and convince yourself that you are at a resonance. In the computer menu, go to “Windows” > “Measure Wave Function”.

Adjust the hemispheres to  $\alpha = 0^\circ$ , and measure the amplitude in steps of  $10^\circ$ . The program converts the angle  $\alpha$  automatically to the polar angle  $\theta$  and plots the absolute of the amplitude in a polar plot. Use the function “complete by symmetry” to complete the figure.


Create polar-plots for the prominent peaks and identify the quantum numbers.

**Analyze data:**


Compare the polar plots you have generated with polar-plots of the Legendre polynomials. Some of them are given below, the others you can visualize with the program PlotYlm.exe.

In case of overlapping peaks, you will find distorted figures, since there are contributions to the wave functions from two different eigenstates with different quantum numbers and symmetries.


**Fig. 2.3:** Plots of the spherical harmonics:

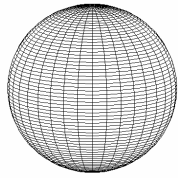


$$Y_0^0(\theta, \varphi)$$

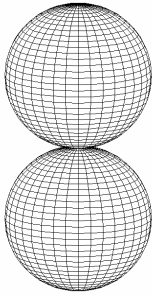


$$Y_1^0(\theta, \varphi)$$

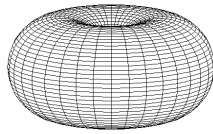


$$Y_1^1(\theta, \varphi)$$

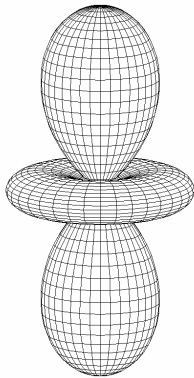


$$Y_2^0(\theta, \varphi)$$

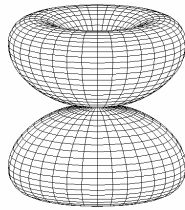


$$Y_2^1(\theta, \varphi)$$

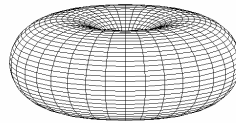


$$Y_2^2(\theta, \varphi)$$

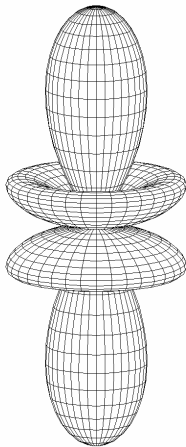


$$Y_3^0(\theta, \varphi)$$

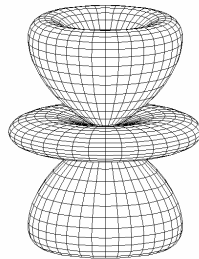


$$Y_3^1(\theta, \varphi)$$

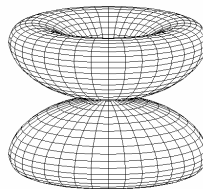

$$Y_3^2(\theta, \varphi)$$

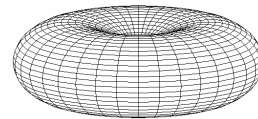



$$Y_3^3(\theta, \varphi)$$

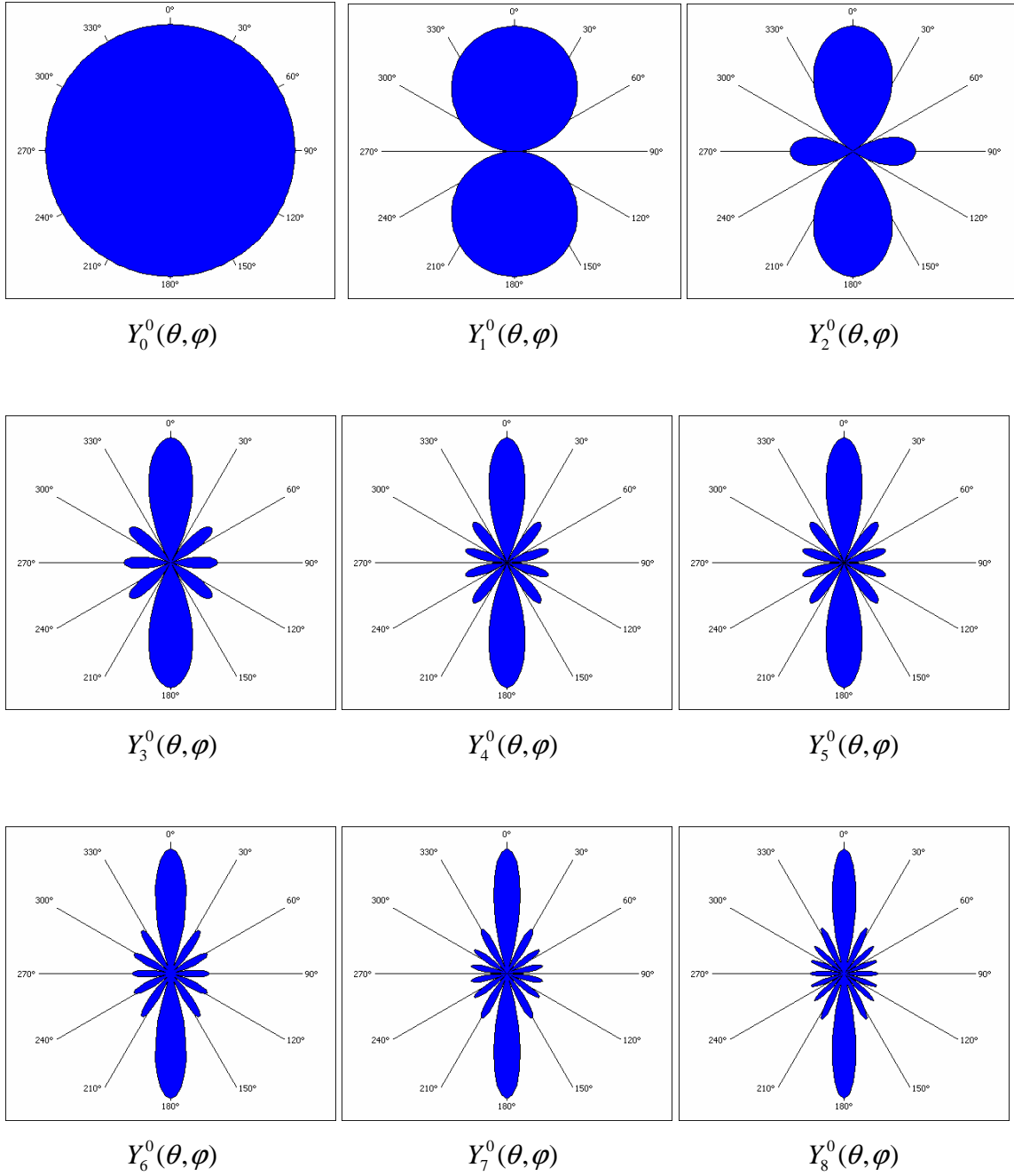**Fig. 2.4:** Cut through the spherical harmonics with magnetic quantum number  $m = 0$ .



**Advisor Information: THE REST OF THIS MANUAL IS FOR THE ADVISOR ONLY**

In the following figures, the wave functions are plotted as measured with the computer. Due to tiny imperfections, the degeneracy of the states with different magnetic quantum numbers is not perfect. Though the lifting of degeneracy is small compared to the width of the peak, it results in modifications of the measured shape of the wavefunction that depend on the exact frequency that has been chosen. Students may want to try measuring at the center of a peak and on its right and left shoulders.




In the following series of figures, the polar plots for most of the peaks below 9 kHz are shown. On the left, the experiment is displayed. On the right, the corresponding calculated Legendre polynomial is shown.

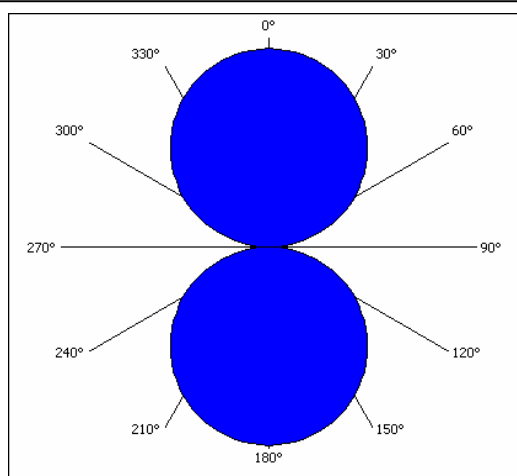
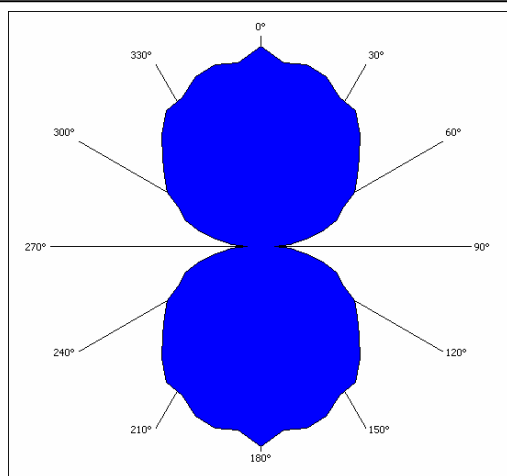


Wavefunction of the peak near 2300 Hz (left) and calculated spherical harmonics for  $l = 1, m = 0$  (right). (Data measured on the right shoulder of the peak).

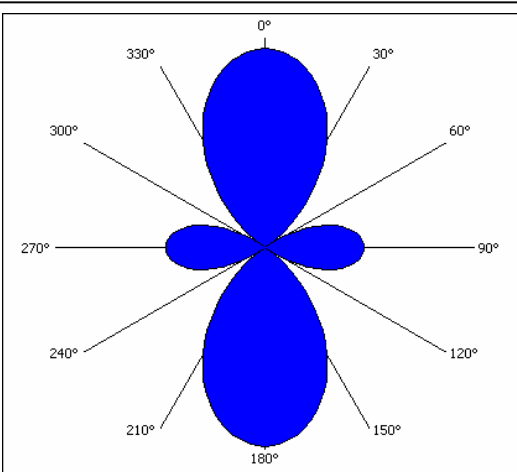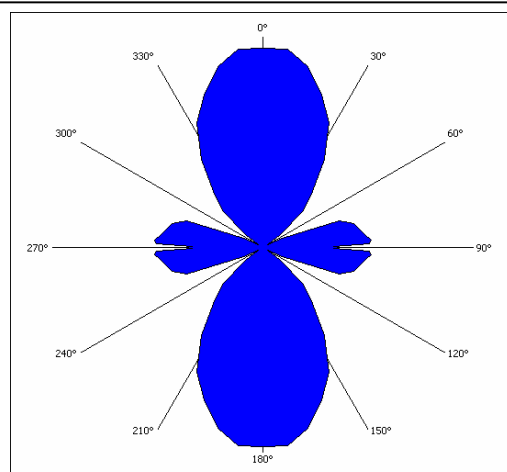



Wavefunction of the peak near 3700Hz (left) and calculated spherical harmonics for  $l = 2, m = 0$  (right). (Data measured on the right shoulder of the peak).

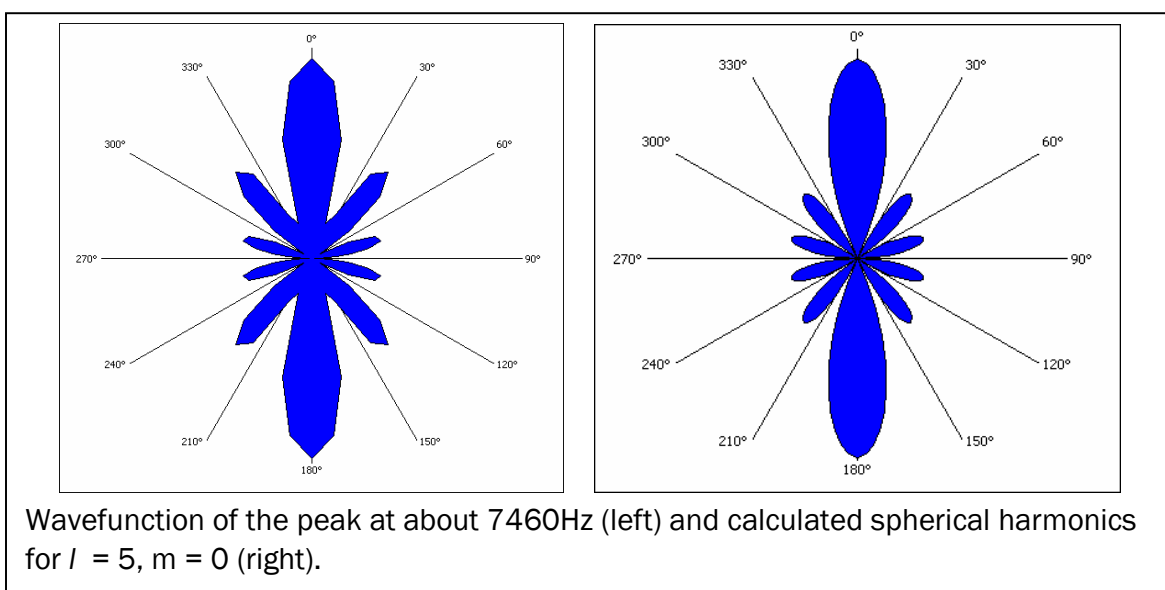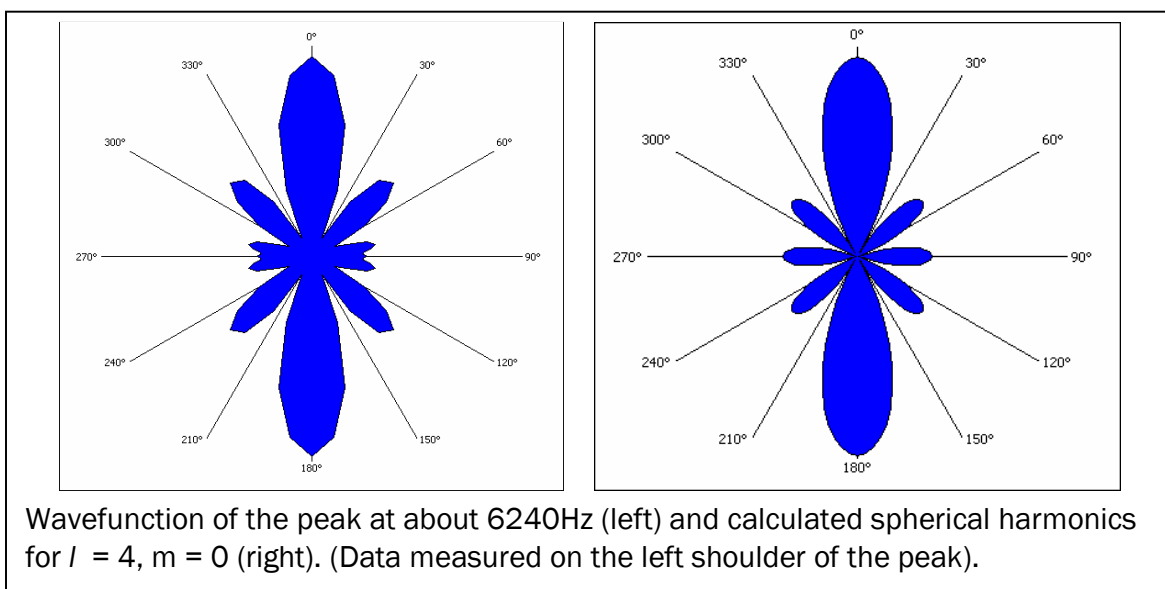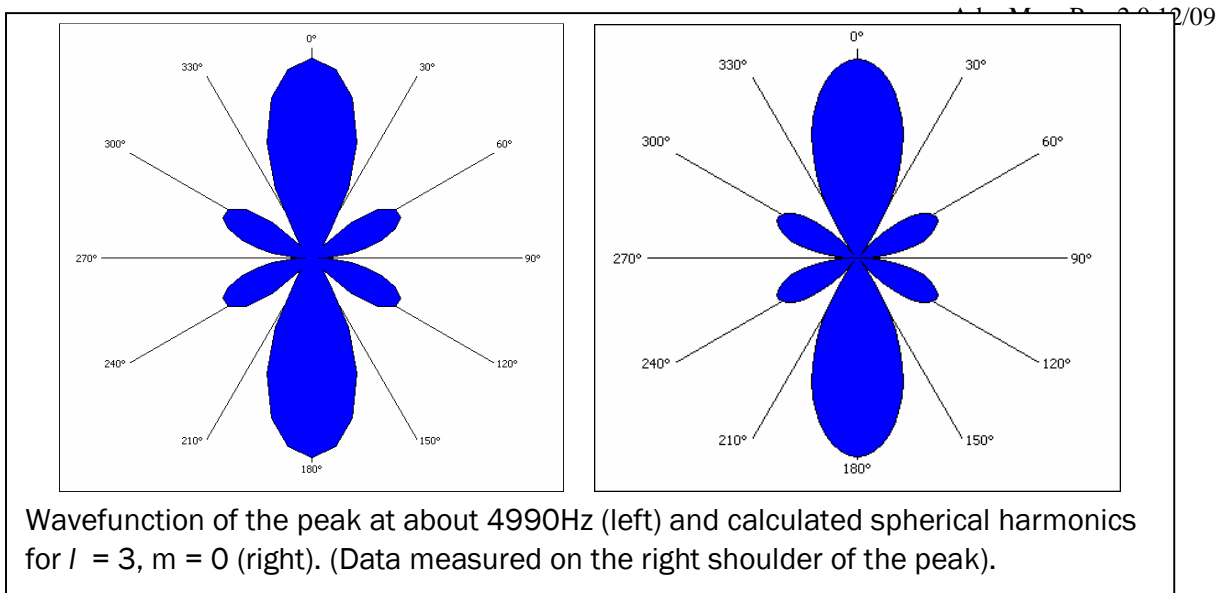

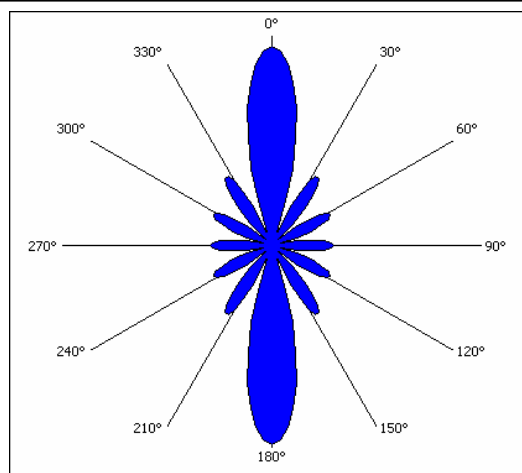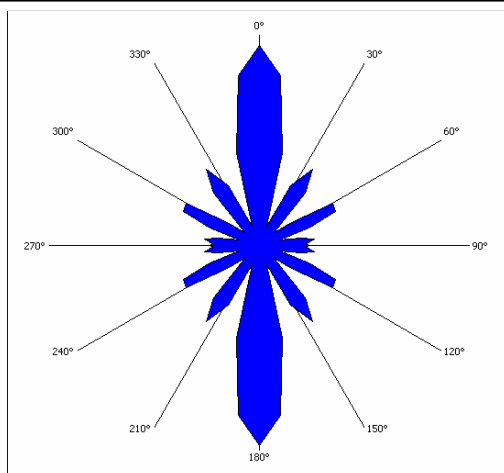


Wavefunction of the peak at about 8650Hz (left) and calculated spherical harmonics for  $l = 6, m = 0$  (right).

Beside the peaks shown above that all belong to the radial quantum number  $n' = 0$  there are peaks belonging to  $n' = 1$ . The wavefunctions of these peaks are shown below.

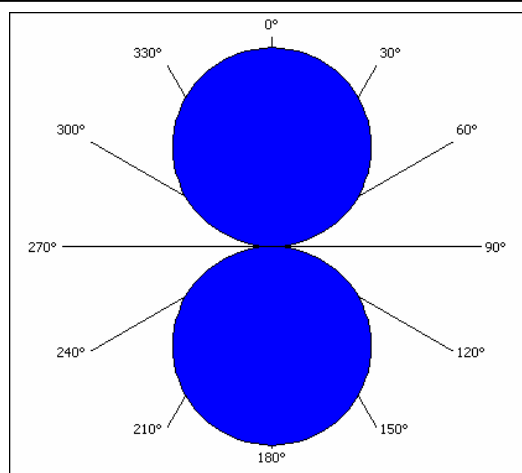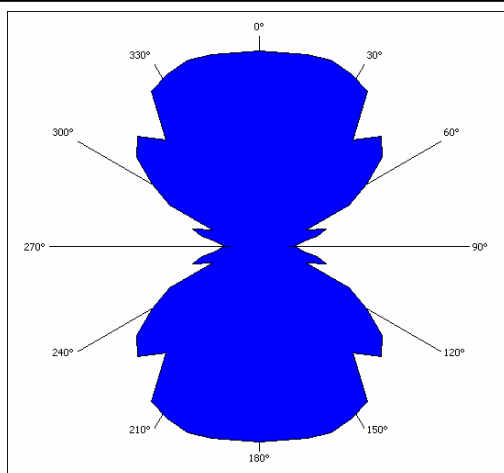


Wavefunction of the peak at about 6570Hz that belongs to the radial quantum number  $n' = 1$  (left) and calculated spherical harmonics for  $l = 1, m = 0$  (right).

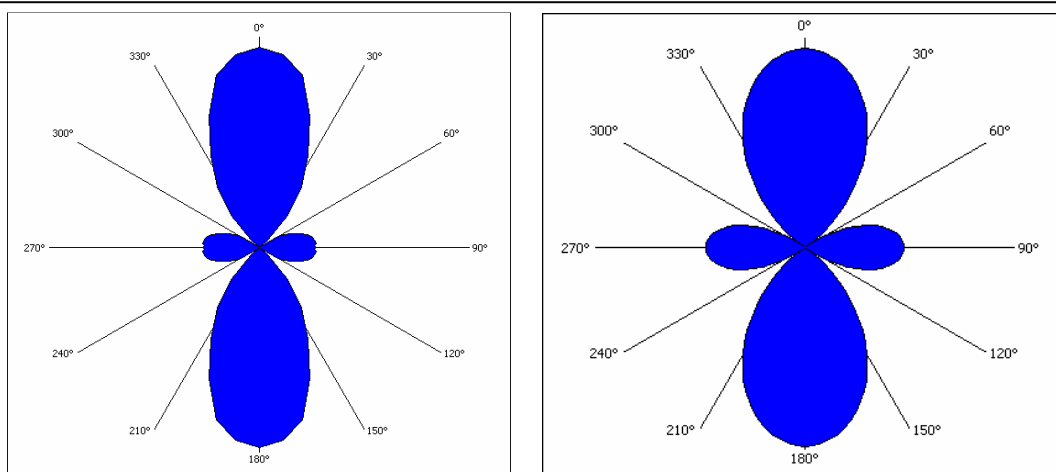


Wavefunction of the peak at about 8040Hz that belongs to the radial quantum number  $n' = 1$  (left) and calculated spherical harmonics for  $l = 2, m = 0$  (right).

After the quantum numbers have been identified by the shape of the wave functions, the peaks can be labeled in the spectrum.

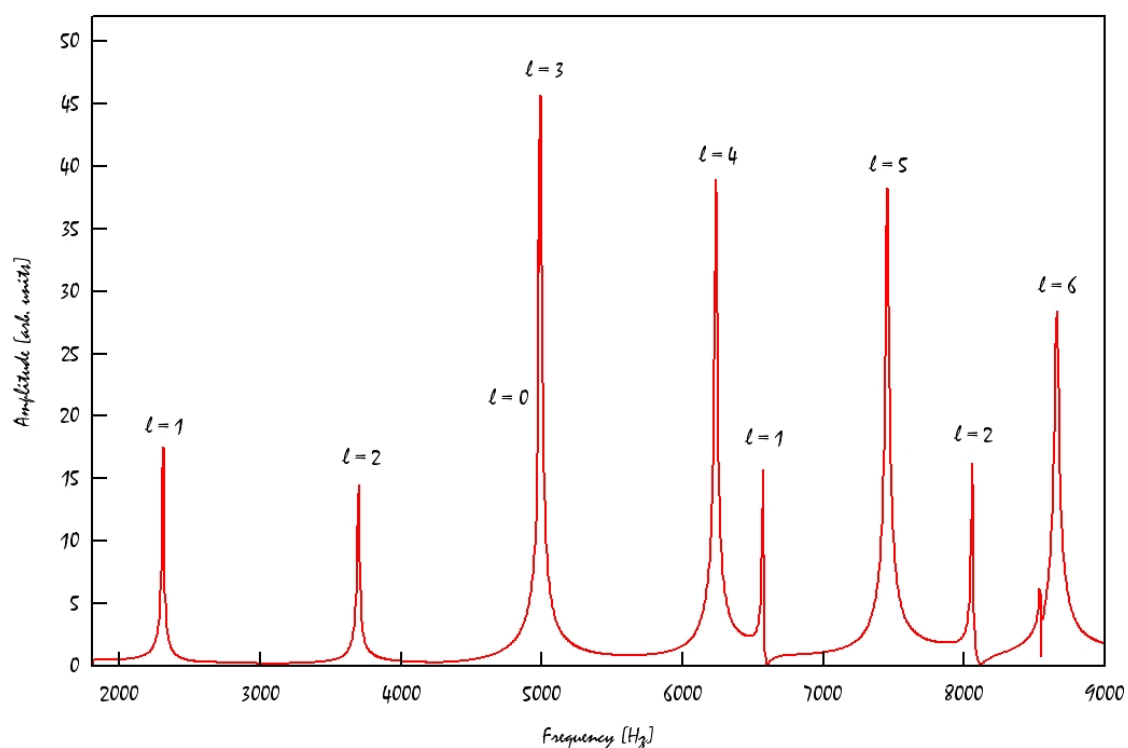

Spectrum measured in the spherical resonator with the peaks labeled by the angular quantum number  $l$ .

Alternatively the peaks can be labeled using the principle quantum number  $n = l + 1 + n'$  as it is convention for the hydrogen atom

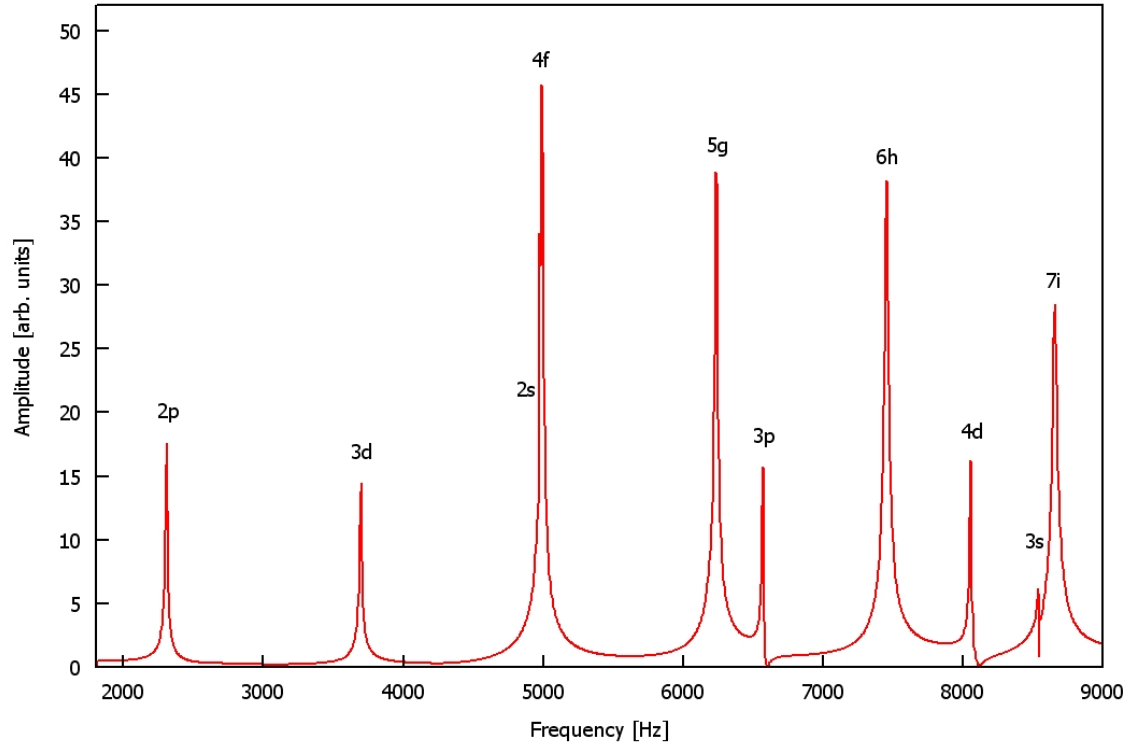

Spectrum measured in the spherical resonator with the peaks labeled in the convention used for the hydrogen atom (principal quantum number  $n$  and a letter to indicate the angular quantum number  $l$  ).